

New !!

New !!

New !!

New !!

EASY CAP EEG-Recording Caps

are well known for

- Superior Signal Quality from Ag/AgCl Sintered Electrodes
- Lowest Impedances together with highest Subject Comfort
- Customer-specific Electrode Arrangements

Serving prime research labs now for 8 years we compiled a

Ready-to-use EEG Recording Cap Set for 21 Channels

incorporating the electrode positions of the

International 10/20 - System.

This package contains everything you need
to start recording right out of the box:

- Cap with 23 electrode positions (10/20-system)
- 26 Ag/AgCl sintered electrodes, number labeled 1 – 26
- 2 Ear-lobe electrodes
- Chest belt set
- All necessary consumables

and as a **totally new feature** to shorten mounting time

- **a quick-connect interface!**

To get all this just order

EC20-#

(Replace the # by the preferred cap size)



EASY CAP EC20

Standardized Ready-to-use EEG Recording Cap Set for 19/21 Channels

consisting of

- | | | |
|----|-----------|-----------------------------------------------------------------------------------------------------------------------------------------------|
| 1 | SU-# | <i>Subtemporal cap for # cm head circumference*
equipped with 23 electrode adaptors A06
according to Montage No. 25 (see graphic)</i> |
| 26 | B10-S-120 | <i>AgCl sintered ring electrode
with 120 cm Heavy-duty leadwire
and 1,5mm touchproof safety socket</i> |
| 1 | C4 | <i>Pair of Nihon-Kohden ear clips
with appropriate electrodes</i> |
| 1 | KV | <i>38-Channel-Electrode board adaptor</i> |
| 1 | KI | <i>38-Channel-Input box</i> |
| 1 | GS | <i>Chest belt set:
1 Chest belt, 2 chin straps, 1 chin belt</i> |
| 1 | V17 | <i>Abralyt 2000 (1 kg)
Chloride-free, abrasive electrolyte gel</i> |

Consumables:

*Alcohol, cotton swabs, syringe,
measuring tape, spare adaptors,
cleaning brush, manual,
channel layout graphic*

The catalogue number of the entire package is EC20-#*.

When re-ordering, please use the catalogue numbers of the single items, e.g.

- | | | |
|---|------|---------------------------------------------------------------------------------------------------------------------------------|
| 1 | SU-# | <i>Subtemporal cap for # cm head circumference*
equipped with 23 electrode adaptors A06
according to Montage No. 25</i> |
|---|------|---------------------------------------------------------------------------------------------------------------------------------|


Our full product list can be seen at www.advancedmedicalequipment.com

* Cap sizes are given in cm head circumference:

- Adult caps: 62, 60, 58, 56, 54
- Children caps: 54, 52, 50
- Infant caps: 36, 38, 40, 42, 44, 46, 48

Montage No. 25

23 Equidistant electrode sites according to Int. 10/20-System
For Reference, FCz or TP9/TP10 or ear lobes could be used
AFz could serve as Ground
Spare electrodes can be placed on skin for use as EOG, ECG, etc.

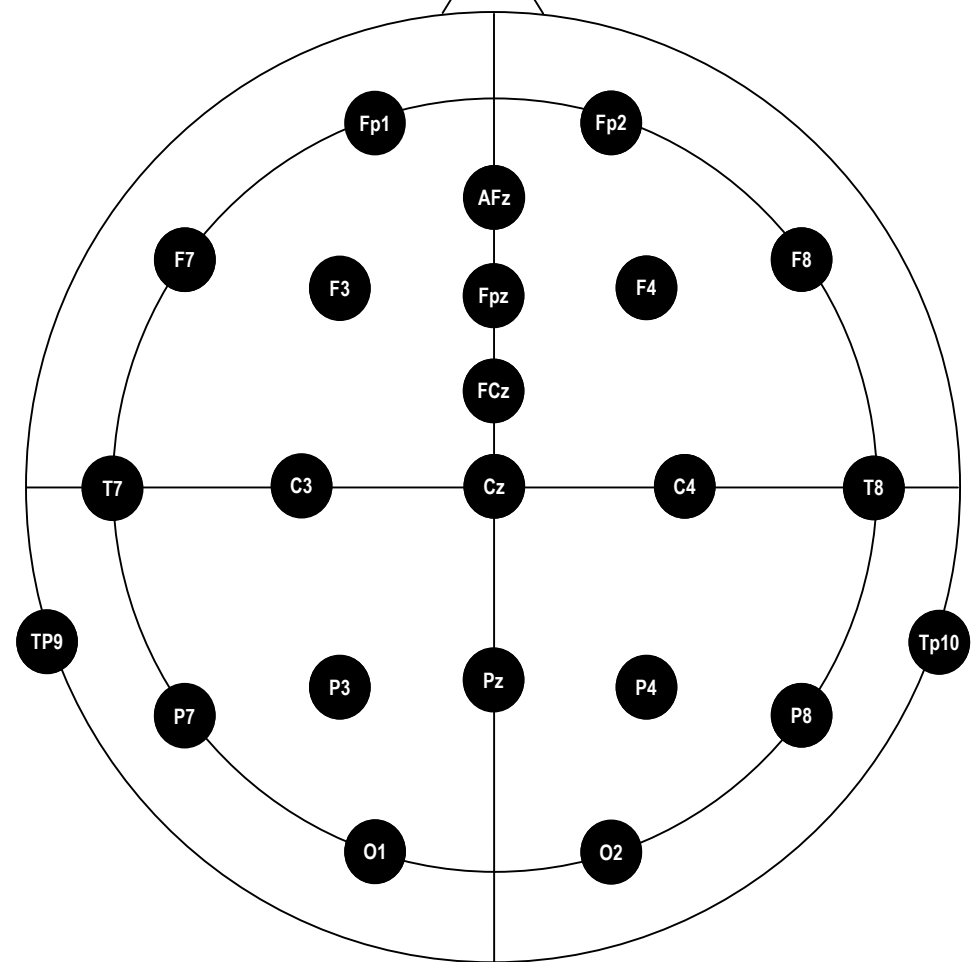

Presented by:

advanced
MEDICAL
equipment ltd

Tel: 01403 260156

Fax: 01403 260175

info@advancedmedicalequipment.com