

New !!

New !!

New !!

New !!

## **EASY CAP EEG-Recording Caps**

are well known for

- Superior Signal Quality from Ag/AgCl Sintered Electrodes
- Lowest Impedances together with highest Subject Comfort
- Customer-specific Electrode Arrangements

Serving prime research labs now for 8 years we compiled a

### **Ready-to-use EEG Recording Cap Set for 32 / 40 Channels**

incorporating the most important electrode positions  
for all EEG machines with 32 or 40 Channels.

This package contains everything you need  
to start recording right out of the box:

- Cap with 36 equidistant electrode positions (10%-system)
- 38 Ag/AgCl sintered electrodes, number labeled 1 – 38
- 2 Ear-lobe electrodes
- Chest belt set
- All necessary consumables

and as a **totally new feature** to shorten mounting time

- **a quick-connect interface!**

To get all this just order

## **EC40-#**

(Replace the # by the preferred cap size)



## EASY CAP EC40

### Standardized Ready-to-use EEG Recording Cap Set for 32/40 Channels

consisting of

- 1 SU-# *Subtemporal cap for # cm head circumference\* equipped with 36 electrode adaptors A06 according to Montage No. 24 (see graphic)*
- 38 B10-S-120 *AgCl sintered ring electrode with 120 cm Heavy-duty leadwire and 1,5mm touchproof safety socket*
- 1 C4 *Pair of Nihon-Kohden ear clips with appropriate electrodes*
- 1 KV *38-Channel-Electrode board adaptor*
- 1 KI *38-Channel-Input box*
- 1 GS *Chest belt set:  
1 Chest belt, 2 chin straps, 1 chin belt*
- 1 V17 *Abralyt 2000 (1 kg)  
Chloride-free, abrasive electrolyte gel*

#### Consumables:

*Alcohol, cotton swabs, syringe,  
measuring tape, spare adaptors,  
cleaning brush, manual,  
channel layout graphic*

The catalogue number of the entire package is EC40-#\*.

When re-ordering, please use the catalogue numbers of the single items, e.g.

- 1 SU-# *Subtemporal cap for # cm head circumference\* equipped with 36 electrode adaptors A06 according to Montage No. 24*


Our full product list can be seen at [www.advancedmedicalequipment.com](http://www.advancedmedicalequipment.com)

\* Cap sizes are given in cm head circumference:

- Adult caps: 62, 60, 58, 56, 54
- Children caps: 54, 52, 50
- Infant caps: 36, 38, 40, 42, 44, 46, 48

#### Montage No. 24

34 Equidistant electrode sites, chosen from 10%-system  
For Reference, FCz or TP9/TP10 or ear lobes could be used  
AFz could serve as Ground  
Spare electrodes can be placed on skin for use as EOG, ECG, etc.

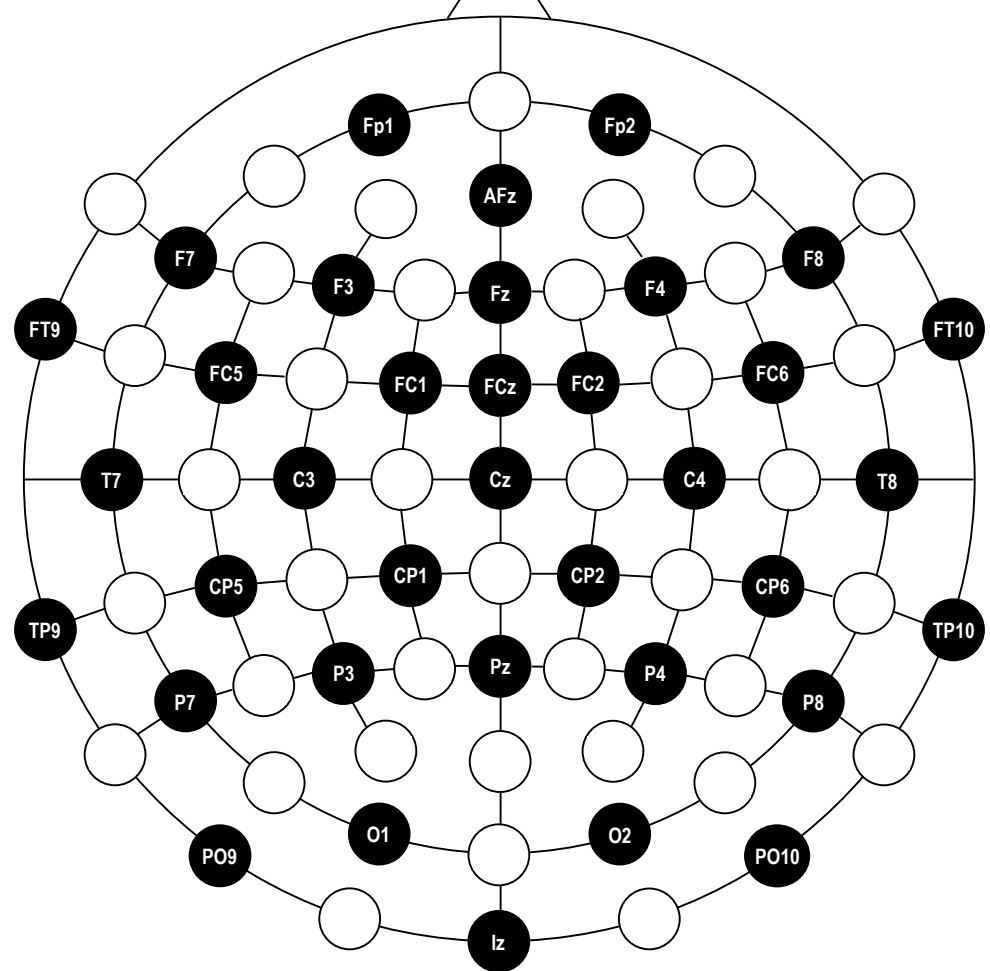

Presented by:

advanced  
**MEDICAL**  
equipment ltd

Tel: 01403 260156  
Fax: 01403 260175  
[info@advancedmedicalequipment.com](mailto:info@advancedmedicalequipment.com)